

televic

AQURA Essentials, the core of nurse call

Our upgradable and versatile
entry level nurse call system,
where quality meets price.



HEALTHCARE

Innovative concept:

- › Very easy to install: no prior nurse call knowledge needed
- › Auto detection: no configuration needed
- › Units: call buttons and outlets
 - Clickable for swift installation
 - Concave design avoiding false alarms
 - No dirt gaps
- › Intelligent fall back to room mode in case of system failure
- › Identification of call per bed
- › Different alarm types:
 - Normal • Pull out • Sanitary • Monitoring
 - Assistance • Sanitary assistance • Medical



Reliable
and safe nurse call system



User-friendly
for nurse, technician and patient/resident



Innovative
thanks to the advanced technology and design



Robust
high quality materials designed for care and cure environments



Upgradable
to the full AQURA care platform




AQURA Essentials

Discover the upgradable and versatile entry level nurse call system of Televic Healthcare. A budget friendly answer to the current challenges of healthcare organizations, with advantages for every stakeholder.

Choosing for AQURA Essentials means:

- › Plug-&-play nurse call system
- › No server needed
- › Capacity up to 40 rooms per ward controller
- › Alarm distinction per bed
- › Fully customizable signage for alarms
- › Nurse call forwarding within your ward
- › Intuitive touch nurse station
- › Hands-free digital intercom
- › UV resistant and disinfectant

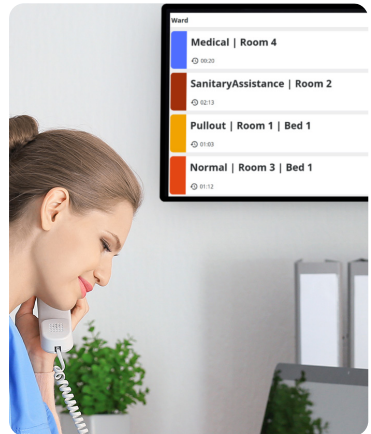

Protect your investment and upgrade to the AQURA Care platform

Over time the needs of your organization change, AQURA grows with you. At any time you can upgrade to the AQURA care platform by simply adding a server and the licenses needed. No change of hardware, nor cabling changes are needed:

- › Mobile alarming on smartphones
- › Intelligent alarm distribution
- › Insights through dashboarding
- › Integration with third party applications & sensors



Smart intercom unit

- › Full duplex digital audio
- › Crystal clear: echo cancellation & noise suppression



Smart single outlet with call and cancellation button

- › Money saving magnetic connector with spring contacts
- › Easy bed assignment
- › Room- and reading light support
- › Ability to detect a dry contact alarm
- › Blinds control in combination with IO extender




Push pears & monitoring cable

- › 1-button: nurse call only, 5-button: nurse call with room- and reading light control
- › Finding light to find the push pear in the dark
- › Flame-retardant and anti-static
- › IP67 rating



Nurses can focus on the residents or patients, thanks to the intuitive user interface.

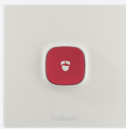


Call unit for normal, sanitary or assistance calls

- › Integrated finding light (LED)
- › Nurse symbol for easy differentiation
- › Reassurance light




Components of AQURA Essentials

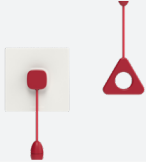

Pull cord unit for sanitary, normal or assistance calls

- › Integrated finding light (LED)
- › Nurse symbol for easy differentiation
- › Reassurance light
- › Optionally: IP54

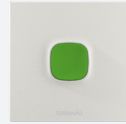

Nurse presence and call cancellation unit

- › Integrated finding light (LED)
- › Nurse call forwarding by default (buzzer)




Overdoor light with integrated room hub

- › Clear distinction of alarms
- › LED technology
- › 5 colours:
 - Red (first row)
 - Green (second row)
 - Yellow or white (third row)
 - Blue or white (fourth row)
- › Intelligent fallback to local room mode in case of system failure



An innovative system

AQURA Essentials allows you to optimize your investment over time: it is a scalable solution that grows with you as the needs of your organization change.

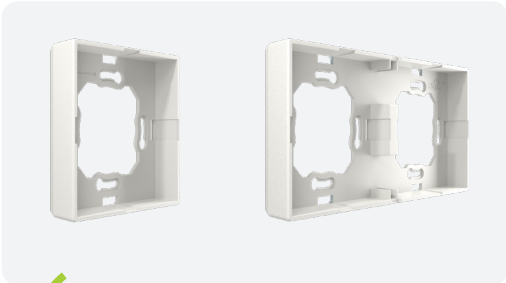



Smart unit for medical alarm

- › Blue button with integrated finding light for medical or reanimation alarm
- › Yellow button to reset the alarm
- › Integrated finding light (LED)
- › Reassurance light



Signage

The AQURA Essentials web based application visualizes the 4 highest active alarms, in real time, on the display of your choice:

- › Show the duration each call has been active
- › Repeat a sound rhythm of the call with the highest priority
- › Mute the sound rhythm in case it is disturbing



Nurse station touch phone

A user-friendly phone with an interactive Android touch interface, showing all active alarms of the ward and enabling intercom conversations with the patient room or at bed level.

- › One device for intercom, interactive alarm monitoring and handling
- › Easy setup to local language for alarm types and user interface
- › Automatic prioritization of alarms
- › Click on the alarm and directly setup an intercom call
- › Waiting time indication
- › Day and night regime for audio volumes
- › 24-hours history of calls
- › Runs on Android, ensuring a smooth and intuitive handling
- › Available in two sizes: 7 or 8 inch



Smart ward controller

The ward controller is the control and communication unit:

- › Capacity of up to 40 rooms per ward
- › Easy room assignment: no programming needed
- › Fallback mechanism to room mode in case the ward controller gets isolated from the system
- › Ability to split the ward in 2 branches for extra reliability
- › DIN rail mountable



The modular and clickable units can easily be installed and are plug-and-play.



Receive all events and alarms on your smartphone by upgrading to the AQURA Care platform



IO extender

IO extender for blind control and increased cabling flexibility when using light-or-blind control.

televic



The safety and happiness of your residents, patients and nurses is our main priority.



HEALTHCARE



Televic Healthcare
Leo Bekaertlaan 1
B-8870 Izegem, Belgium
+32 51 30 30 45
global.healthcare@televic.com
www.aqura-essentials.com